JAY®

LitePro™ CUSHION

JAY® LitePro™ Cushion

Important Consumer Information

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NOTE:

- Check all parts from shipping damage. In case of damage DO NOT use. Contact Carrier/Sunrise Medical for further instructions.
- Sunrise Medical recommends that a clinician such as a doctor or therapist experienced in seating and positioning be consulted to determine if a JAY LitePro Cushion is appropriate.
- Cushions should only be installed by an authorized Sunrise Medical supplier.

⚠ GENERAL WARNINGS

- DO NOT install this equipment without first reading and understanding this manual. If you are unable to understand the Warnings, Cautions and Instructions, contact a qualified clinician or supplier - otherwise injury or damage may occur.
- Prior to prolonged sitting, any cushion should be tried for a few hours at a time while a clinician inspects your skin to ensure that red pressure spots are not developing. You should regularly check for skin redness. The clinical indicator for tissue breakdown is skin redness. If your skin develops redness, discontinue the use of the cushion immediately and see your doctor or therapist.
- All JAY Cushions are designed to help reduce pressure. However, no cushion can completely eliminate sitting pressure or prevent pressure sores.
- JAY Cushions are not a substitute for good skin care including; proper diet, cleanliness, and regular pressure reliefs.
- Installing a cushion on a wheelchair may affect the center of gravity of the wheelchair and may cause the wheelchair to tip backwards, potentially resulting in injury.
- Always assess for the potential need for antitippers or amputee axle adapter brackets to be added to the wheelchair to help increase stability.
- Hook and Loop material should not come into prolonged contact with moisture as this may degrade the adhesive and lead to a failure of the Hook and Loop material.
- Allow the cushion to warm or cool to room temperature before using. Sitting on extremely hot or cold surfaces can cause skin damage.
- OBSTRUCTIONS: DO NOT place any obstructions between the user and the cushion because it will reduce product effectiveness.
- The pods contained in the insert are determined by the clinician during a seating evaluation. They should not be adjusted without a clinician directive as they can alter the pressure relieving properties and lead to injury.
- CUSHION AND COVER ORIENTATION: Product must be used with the insert facing up. If the cover is not used correctly it may reduce or eliminate the cushion's benefits and could increase risk to the skin and soft tissue.
- Avoid sharp objects or exposure to excessive heat (above 108°F / 42°C) or open flame.
- Freezing temperatures (less than 32° F / 0° C) can make the insert more firm and may affect its pressure relieving qualities. When using the cushion in freezing temperatures, perform skin checks and pressure reliefs more often to ensure proper skin protection.
- Exposure to heat and repeated pressure may make the insert deform slowly over time, compromising its pressure relieving capabilities. Make sure to check your cushion monthly, or more often in warmer conditions, for signs of deformation.

Word	Definition
 DANGER!	Advice to the user of potential risk of serious injury or death if the advice is not followed
 WARNING!	Advice to the user of a risk of injury if the advice is not followed
 CAUTION!	Advice to user that potential damage to equipment may occur if the advice is not followed
NOTE:	General advice or best practice
	Reference To Additional Documentation

INTENDED USE**Intended Use**

The JAY® LitePro™ is a wheelchair cushion that aids in the stabilization of the pelvis and torso of the wheelchair user and properly positions the wheelchair user's body while seated.

The JAY LitePro Cushion is designed for clinicians to use with clients who need a high level of pressure distribution, stability and skin protection. The multizone elastomer insert has interchangeable pressure adjustment pods with two levels of firmness to adjust the level of pressure management and accommodate pelvic obliquities. Pods should only be changed at the direction of your clinician.

Maximum user weight for this product is 150 Kg.

The intended lifetime of this product is 5 years.

Before use please consult the attending doctor, as to whether you can/should use the cushion and if so, for how long.

JAY LITEPRO COMPONENTS

FITTING A PRESSURE DISTRIBUTING CUSHION

Select a cushion to match your hip width when in the seated position.

- The cushion length should extend to within 25 mm - 50 mm from the back of the knees.
- Place the cushion to the back of the wheelchair with the material identification tags and insert toward the rear. When fit correctly, your hips should be to the back of the chair. Your ischial tuberosities (sit bones) should be centered on the insert in the pelvic loading area of the cushion.
- Proper footrest adjustment may enhance sitting comfort and help lower peak sitting pressures.
- Adjust the footrests so that your legs rest firmly but comfortably on the cushion.

NOTE: Do not use the wheelchair's seat width as a guide in determining the appropriate cushion width. Improper cushion fit may be detrimental to the performance of a pressure distribution cushion and its support of the body.

Avoid using the cushion on overly stretched seat upholstery. Stretched upholstery does not provide adequate support for the cushion and improper positioning may result. Replace stretched upholstery or purchase the JAY Adjustable Solid Seat.

Contact an authorised supplier for a free demonstration of these products.

CUSHION ADJUSTABILITY

The JAY LitePro Cushion offers adjustability to accommodate for pelvic obliquity through interchangeable pressure adjustment pods with differing levels of firmness, which are located on the underside of the insert (Fig. 1.0). The determination of the correct combination of pods should be made by your clinician during evaluation. Over time, as the user adapts, the pods may be removed or added to allow for maximum cushion to user contact. Pods should only be changed in the direction of your clinician.


Pods aligned flush in insert

Fig. 1.0

! WARNING!

Place the insert in the base, making sure it is aligned properly in the grooves of the base (Fig. 2.0). A misalignment or creation of an edge could cause an unwanted pressure point that may lead to injury.

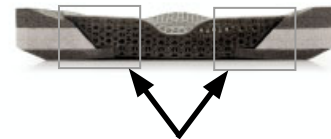Align insert into grooves of base

Fig. 2.0

EASY MAINTENANCE & CLEANING

Maintenance

Monthly cleaning and regular maintenance may help extend the life of your cushion.

Check the cover for tears and excessive wear and replace if ripped, torn, or otherwise not fully functional. Inspect the insert and pods for tears, deformation, or discoloration.

Check the foam base to ensure foam consistency. Excessive wear and tear could compromise the performance of this cushion.

Check the insert and pods for tears, deformation, or discoloration.

If you notice any inconsistencies, please reach out to your local dealer to inspect the components for possible replacement.

To clean cover

Washing:

- Remove the cover from the foam base and turn inside out.
- Machine wash in warm water and drip or tumble dry on low heat (follow instructions from cleaning chart using common laundry detergent).
- No ironing required.

Cleaning

- Remove the cover from the foam base
- Spray the surface with a 70% Alcohol based or 0.1% Chlorine solution
- Let it dry
- Wipe off completely the solution with a moistured, non-abrasive cloth with clear water
- Let it completely airdry

⚠ CAUTION!

- Do not dry clean covers or use industrial washers and dryers to clean cover.
- Do not steam autoclave.
- Do not bleach.

To refit cover after washing/cleaning:

- Reinstall the cover on the foam base.
- Ensure the back of the cover matches up to the back of the foam base.

NOTE: If utilized in an industrial setting, write the resident's name with a permanent marker on the content label of the covers. This will help to ensure the return of the covers after washing.

Use only JAY designed cover on JAY cushions. The cover is an important part of allowing immersion into the insert which redistributes the pressure. Changing the cover will alter this design and its efficacy.

To clean insert

- Remove the cover from the foam base.
- Remove insert from the inner cover and handwash with warm water and soap, or place in the dishwasher, without any pressure on it.
- Let completely air dry before reassembly.
- Reinstall the insert on the base and reinstall the cover over the foam base.
- Ensure the back of the cover matches up to the foam base.

NOTE: Normal household cleaners such as detergent, sanitizer, and bleach are all safe to use on the insert. Do not use strong alcohol-based cleaners, such as hydrogen peroxide, to clean the insert. Avoid exposing the insert to excessive oils such as lotions. Exposure to excessive amounts of strong acids and oils can make the elastomer material soften, compromising the pressure-relieving abilities of the insert.

NOTE: If you do use bleach on the insert, make sure it is completely rinsed and dried before putting back together with the other components of the cushion as bleach will degrade the foam base and cover of this cushion.

Do not wash a damaged insert.

To clean foam base

- Remove the cover and the insert from the foam base. Wipe lightly with a damp cloth.
- Wipe off the base with a clean cloth. Let completely air dry before reassembly.
- Reinstall the 3D insert on the base and reinstall the cover over the foam base.
- Ensure the back of the cover matches up to the back of the foam base.

⚠ CAUTION!

- Do not use soap.

If at any time, the cover, foam base, insert, or any other component appears to be wearing or you notice a degradation in the texture, contact your supplier/Sunrise Medical for evaluation and further instructions.

NOTE: Consult the "Approved Cleaning Procedures for JAY Cushions" chart.

Fig. 3.0		Laundry Care Symbols	
	Machine wash		Do not iron
	Only non-chlorine bleach, when needed		Do not dry clean
	Drip dry		

Reassembling the JAY LitePro cushion

Follow Maintenance & Cleaning instructions as described and then reassemble as follows.

- Place the insert in the foam base, lining up the contours to create a smooth transition between the insert and foam base.

! CAUTION!

FAILURE TO LINE UP THE INSERT AND BASE ACCURATELY COULD RESULT IN PRESSURE INJURIES. (Fig. 4.0, Fig. 5.0)

- Insert the foam base into the cover.
- Ensure that the zipper is at the rear of the foam base where the pelvic loading area is located.

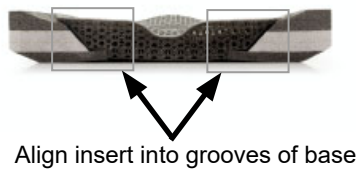
Fig. 4.0

Pods aligned flush in insert

Fig. 5.0

CONDITIONS FOR SAFE STORAGE

- Products should be stored in a dry place: at 20 – 75% Relative Humidity.
- Protect from direct sunlight and dust for long periods.
- Away from direct heat and at a temperature between 5°C to 40°C.
- Position and stack the packaged products according to the «this side up» arrows on the packaging.
- If a package is high and narrow and therefore potentially unstable, ensure it is well secured on the pallet or in the racking, to prevent it from falling down.
- Do not stack anything on top of the packaged products when the «do not stack» mark on the packaging is shown.
- Product should never be stored outside and subject to the elements.

HYGIENE MEASURES WHEN BEING RE-USED

This JAY cushion is designed to be re-used. The frequency of re-use depends largely on the way the cushion has been used previously. It is essential that the cushion is properly cleaned before re-use (see section “Easy Maintenance & Cleaning”). It is strongly recommended to replace the cover when reusing the cushion.

Approved Cleaning Procedures for JAY Cushions										
Product	JAY Lite	JAY LitePro	JAY J3	J2 - J2 DC - GS	JAY Easy Visco	JAY Soft Combi P	JAY Xtreme Active	JAY Basic	JAY Easy	JAY Balance
Outer Cover (Microclimatic and Stretch)										
(1). Wet Wash 70°C 										
Outer Cover (Incontinent)										
(2). Wet Wash 60°C 										
Inner Cover										
(3). Spray and Wipe disinfection or Wet Wash 60°C  	NA	NA	NA	NA	NA	NA		NA	NA	
Cushion core (foam base)										
(4). Steam autoclave up to 105 °C (1) 										
Spray and wipe disinfection 										
(5). Immersion disinfection 										
Cushion elements										
Air element	NA	NA		NA	NA	NA	NA	NA	NA	
Multizone Elastomer insert	NA		NA	NA	NA	NA	NA	NA	NA	NA
Positioning elements	NA	NA			NA	NA		NA	NA	
(Cryo) Fluid pad	NA	NA			NA	NA		NA		
70°C is the maximum recommended temperature. 60°C is the maximum recommended temperature. - The Xtreme Active and Balance inner covers can be washed at 60°C as well as spray and wipe disinfected. Possible discoloration if washed at 60°C but the functional aspect of the covers is not diminished. - An autoclave is a self-locking machine that sterilizes with steam under pressure. - At low temperature (< 30°C) - NA (Not Applicable) as the cushion model does not include the cushion element.										
The product warranty shall not apply to products damaged due to the failure to follow these cleaning instructions or the use of inappropriate cleaning compounds.										
Recommended Disinfectants										
Solutions containing: Alcohol 70% solution base, Chlorine 0,1% solution (must be rinsed always with clear water). Refer to the cleaning instructions procedures included in this manual.										

TECHNICAL DATA FOR THE JAY LITEPRO CUSHION

Model	JAY LitePro
Article Number*	JLPwwddSMEU
Weight in ready to use state (Kg)**	0.92
Anti Slip Coating	Yes
Outer Cover Material	Polyester
3D Printed Insert	EPU 43 Elastomer
Inner Cover Material	N/A
Foam Base Material	Polyethylene (PE) Foam
Foam Density	32 and 64 kg/m ³
Compression Hardness	N/A
User Min Weight (Kg)	No minimum weight
User Max Weight (Kg)	150
Spatial Dimension in ready to use state*	Seat Width range (340-500 mm), Seat Depth Range (360-500 mm)

**For the complete offered sizes (Seat Widths - Seat Depths) and different components, please refer to the Order Form of the product. The scope of product delivery will vary according to the clinical assessment and prescription made by a clinician.*


***Data based on an average configuration*

N/A: (Not Applicable) as the cushion model does not include the element


WARRANTY

JAY LITEPRO CUSHION WARRANTY

Each JAY cushion is carefully inspected and tested to provide peak performance. Every JAY cushion is guaranteed to be free from defects in materials and workmanship for a period of 24 months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within 24 months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears or burns, nor to the cushion's removable cover. The removable cover is guaranteed to be free from defects in materials and workmanship for a period of six (6) months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within six (6) months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge. Claims and repairs should be processed through the nearest authorized supplier. Except for express warranties made herein, all other warranties, including implied warranties of merchantability and warranties of fitness for a particular purpose are excluded. There is no implied warranty beyond what is contained herein. Remedies for breach of express warranties herein are limited to repair or replacement of the goods. In no event shall damages for breach of any warranty include any consequential damages or exceed the cost of non-conforming goods sold.

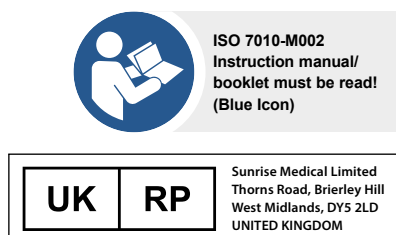Type:	Product Name/ SKU Number	XXXX-XX-XX	Date of Manufacture
LOT	LOT: defined by Julian Date	Part Code	Assembly Part Number
XXX mm	Seat width	MD	This symbol means Medical Device
XXX mm	Seat depth		Manufacturer's address
XXX kg	Max user weight	EC REP	European Authorised Representative
CE	CE Mark	UK RP	UK Responsible Person
UK CA	UKCA Mark	CH REP	Swiss Representative's address
	Consult instructions for use		Importer's address

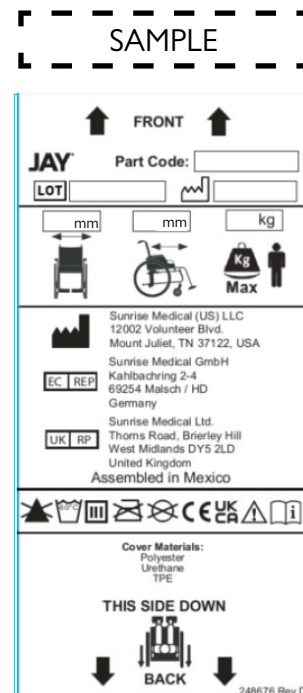The management system of SUNRISE MEDICAL is certified to EN ISO 13485 and ISO 14001.

We at SUNRISE MEDICAL have been awarded the ISO-13485 certificate, which affirms the quality of our products at every stage, from R & D to production. This product meets the requirements in accordance with EU and UK regulations. Options or accessories shown are available at extra cost.

Indications

The varieties of cushions fitting variants, as well as the different positioning components to deliver higher positioning and skin protection benefits according to each user requirement and indication, mean that it can be used by a range of users with different positioning and skin protection needs and indications as follows [R7-R14]:

- Spinal Cord Injury - MS
- Neurological disease pattern (hemiplegia, M.Parkinson)
- Extremity defect deformity - Joint contracture
- Amputation - Tetraplegia

Contraindications

The cushion for a wheelchair shall not be used in case of [R3, R4]:

- Redness developed in the skin contact area
- Bottoming out present in the JAY Fluid Insert
- Decubitus

NOTE: General user advice. Not following these instructions may result in physical injury, damage to the product or damage to the environment!

Notice to the user and/or patient: any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

As the manufacturer, SUNRISE MEDICAL, declares that this product conforms to the Medical Device Regulation (2017/745).

As the manufacturer, SUNRISE MEDICAL, declares that the product conforms to the UK Medical Devices Regulation 2002 No. 618.

B4Me special adaptations

Sunrise Medical strongly recommends that in order to ensure that your B4Me product operates, and performs as intended by the manufacturer; all the user information supplied with your B4Me product is read and understood, before the product is first used. Sunrise Medical also recommends that the user information is not discarded after reading it, but it is kept safely stored for future reference.

Medical Device Combinations

It may be possible to combine this Medical device with one or more other Medical Device or other product. Information on which combinations are possible can be found at www.SunriseMedical.co.uk. All combinations listed have been validated to meet the General Safety and Performance Requirements, section 14.1 of the Medical Device Regulation 2017/745. Guidance on the combination, such as mounting, can be found at www.SunriseMedical.co.uk.

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©2026 Sunrise Medical (US) LLC
Sunrise Medical B.V.
Sunrise Medical HCM B.V.
Sunrise Medical S.A.S