

# JAY Zip Cushion

A comfortable, clinically effective cushion designed uniquely for kids.



The JAY Zip cushion was designed through the use of anthropometric data so they are specially sized to meet a pediatric client's unique needs. It is a dual layered contoured foam base that provides high skin protection and mild positioning in a comfortable, effective and convenient cushion.

## Comfortable

- Dual layered contoured foam base
- 3DX™ fabric spacer in the outer cover promotes air movement
- Soft and stretchable outer cover (also available in pink and blue)



## Convenient & Antimicrobial

### Outer Cover

- Machine washable and quick to dry
- Silver impregnated X-static® stretch fabric naturally inhibits bacterial growth
- Extra layer of microclimatic spacer fabric for comfort
- Safe non-toxic material



**X-static**  
SILVER SYSTEMS™



Outer and inner cover

### Inner Cover

- Easy to wipe off
- Aquaguard™ zipper resists moisture
- Anti-wicking seam thread protects foam base



Dual layered contoured foam base

Easy washing and cleaning


Cushion Width													
Cushion Depth	cm	20	22	24	26	28	30	32	34	36	38	40	
	20												
	22												
	24												
	26												
	28												
	30												
	32												
	34												
	36												
	38												
	40												
	42												
	44												
	46												

Note: This information is for reference only. All prescriptions should be based on clinical assessment findings.

Height (in cm)		Maximum Height
Front	Rear	
5	4,5	6

weight cushion 30x30cm	0,9 kg
max. user weight	75 kg
Guarantee	2 years

These sizes are available, please see your order form for modifications or special sizes.

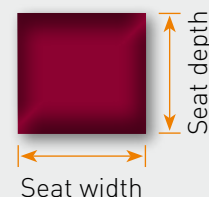

### Sunrise Medical GmbH

Kahlbachring 2-4  
69254 Malsch/HD  
Tel.: +49 (0) 72 53 / 9 80-0  
Fax: +49 (0) 72 53 / 9 80-222  
www.SunriseMedical.eu