



JAY®

LITE PRO

SUPERIOR CLINICAL SEATING – NOW AND THEN

Founded in 1982, JAY revolutionised wheelchair seating with a mission to

- Protect skin
- Support posture
- Enhance mobility

JAY delivers superior clinical seating through developing materials, science, and innovation.

This mission still drives us today. Highly active self-propellers need seating that is more supportive and boosts propulsion. A lightweight cushion that provides reliable support during every transfer, every push, every challenge.

Empower active users and athletes to pursue and achieve their goals with the **JAY LitePro cushion**.

The JAY LitePro cushion provides uncompromising performance in an ultra-light, breathable, and low-profile design.

At just 0.9 kg (2 lbs) and 51 mm (2"), enjoy adjustable skin protection, exceptional stability, and effective pressure management without compromise.





SUPERIOR CLINICAL SEATING – NOW AND THEN

LIGHTWEIGHT AND LOW PROFILE

At 0.9 kg (2 lbs) and only 51 mm (2"), the JAY LitePro is **lightweight and ultra low profile** for improved transfer, agility and performance.

ADJUSTABLE PRESSURE RELIEF AND STABLE SUPPORT

With an contoured pelvic well and a shaped elastomer insert, the JAY LitePro is desgind for perfomance and optimal stability. With an adjustable level of pressure management, the JAY LitePro combines **pressure management with high stability**.

MICROCLIMATE CONTROL

With a perforated closed-cell foam base and breathable elastomer well insert, the JAY LitePro maximises airflow to **keeping the client cool and dry** for stronger, healthier skin.

DURABLE AND EASY TO USE

Long-lasting performance paired with easy, **low-maintenance** materials makes the JAY LitePro perfect for those on the go.

SUPPORT FOR UNCOMPRIMISING PERFORMANCE

INNOVATIVE MULTI-ZONE DESIGN

Innovative elastomer insert is built with different zones.

Each individual zone is engineered with its own density profile, providing support where the body needs stability and pressure relief to protect bony prominences.

The different density zones use anthropometric data, measurements of the size and proportions of the human body to determine optimal zone dimensions.



Multi-density zones

EXCEPTIONAL STABILITY IN A LOW PROFILE CUSHION

The LitePro offers exceptional stability in a low-profile cushion without compromising pressure relief.

Contoured, multi-layer foam delivers postural support and stability, while the immersive insert balances pressure relief with targeted support.

Together, they promote stable pelvic positioning to **improve propulsion, posture and transfers.**



Contoured, multi-layer foam



Highly breathable for less moisture build up improving microclimate and lowering the risk of skin breakdown



Ultra-lightweight cushion at under 1kg (2 lbs) light to push, easy to handle and transport



Simple to clean, low maintenance materials for long-lasting performance



Immersive zones for targeted pressure redistribution, reducing risk at bony areas




Contoured, multi-layer foam for postural support and stability



Contouring and well shaping for high stability in a low-profile cushion



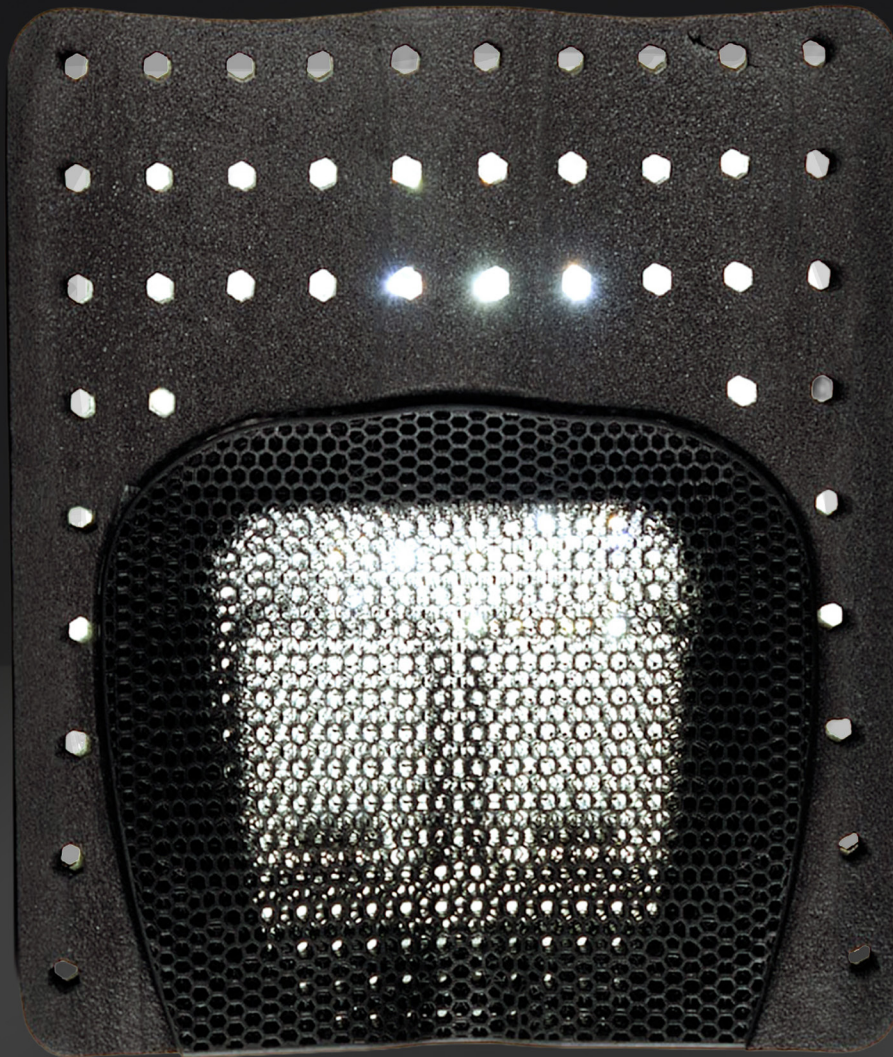

KEEPS SKIN DRY AND COOL


360° BREATHABLE

The elastomer insert paired with the perforated foam base provides continuous airflow throughout the cushion, helping keep temperatures lower where it matters most.

Horizontal and vertical channels integrated into the perforated foam base enhance cooling and moisture evaporation, maximising airflow.

Air Exchange cover promotes additional ventilation at the skin interface, helping keep the user cooler and drier, promoting healthier skin.

The LitePro keeps the most active users cool and dry, providing enhanced comfort and skin protection from best-in-class breathability.

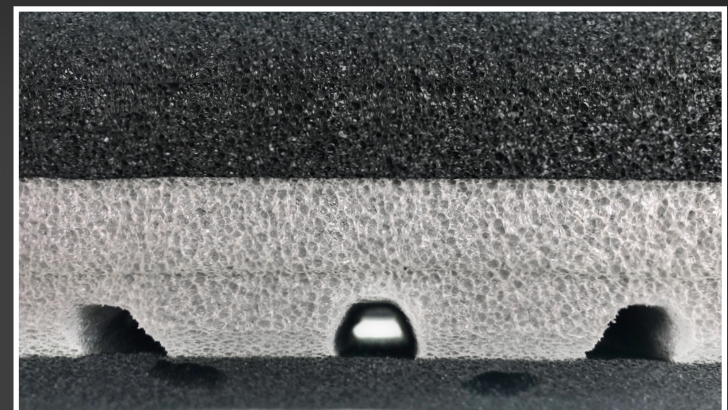

ADJUSTABLE PRESSURE RELIEF

DEEP IMMERSION AND LOW PEAK PRESSURES

Select up to 40 mm of immersion through adjustable pressure management pods, **redistributing pressure away from high-risk bony areas**, for moderate to high skin protection.


The elastomer sits firmly in the LitePro's contoured well, maintaining consistent skin protection even during frequent transfers and active daily movement.

The pressure management pods help **tailor the pressure relief** for an individual addressing a mild obliquity.

COMFORT THAT ADAPTS

The elastomer insert flexes in multiple directions to absorb pressure and impact from every angle, **adapting to the client's unique contours**.

The elastomer offers dampening and impact absorption, creating **even pressure distribution** and a soft, cradled feel for lasting comfort and stability.

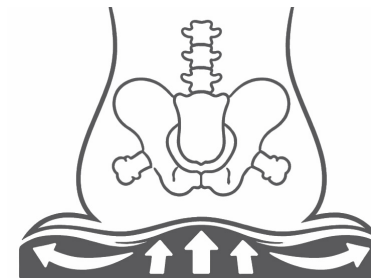

Pressure management pods - available in 2 densities



TECHNICAL DATA	
WEIGHT CAPACITY	150 kg (330 lbs)
WIDTHS	EU 340 mm – 500 mm (13" – 20") US 360 mm – 500 mm (14" – 20")
DEPTHS	360 mm – 500 mm (14" – 20")
CUSHION HEIGHT	50.8 mm (2")
CUSHION WEIGHT	0.90 kg (2 lbs)
COVER MATERIALS	Air Exchange microclimate outer cover with durable side bands and PVC-free non-slip mesh
JAY YOUR WAY	Modifications Available

Sunrise Medical
Thorns Road, Brierley Hill,
West Midlands, DY5 2LD
t: +44 (0) 845 605 6688
e: marketing@sunmed.co.uk
sunrisemedical.co.uk

For further information on the full specification, options, and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.



If you are visually impaired, this document can be viewed in PDF format at sunrisemedical.co.uk

