

QUICKIE[®]
Live without limits.



QS 7^x

/ GO FURTHER. FOLD SMALLER.

QS
SERIES



QUICKIE QS7 X

REDEFINING ULTRA-LIGHTWEIGHT FOLDING PERFORMANCE

“Live without limits” – that is what QUICKIE stands for. The QUICKIE Xenon² was a breakthrough that redefined what a folding wheelchair could be. It set a new benchmark for lightweight design and mobility, helping riders overcome everyday challenges like bulky wheelchairs, complex folding, and heavy lifting.

Today, we build on that legacy with the **QUICKIE QS7 X**. And push it even further.

FOLDS SMALLER. LIFTS LIGHTER. RIDES EASIER. LOOKS BETTER.

Every detail of this next generation of a premium folding wheelchair has been re-engineered to give you one thing: More freedom with less effort.



Scan the QR code to see
the QS7 X Visualizer!





TRANSPORT

EFFORTLESS FOLDING WITH BUILT-IN STABILITY

The QS7 X delivers innovative **FreeFold™** performance for smooth, low-effort folding. Optimized kinematics ensure a controlled, fluid motion for simple, convenient handling.

The QS7 X maintains consistent wheel camber as it folds - **staying stable and upright**. This makes handling and transport easier, especially with frequent folding, so you can place it beside you and keep your hands free - for example, after transferring into your car and before loading it onto the passenger seat.

FOLD-N-LOCK comes standard on the QS7 X, keeping the chair securely folded. With a simple click, it stays compact, with no need to hold it closed during lifting or storage.





UP TO 45% SMALLER THAN STANDARD FOLDING

It slips into tight spaces with ease. Combined with its ultra low weight, the QS7 X makes transport faster, easier, and more effortless than ever.



FITS SPACES WHERE MOST OTHERS DON'T

We pioneered it decades ago. Now we are redefining it: The new **Front Fold** mechanism on the QUICKIE QS7 X folds the front frame, caster link, fork, and caster in one swift movement, cutting transport size dramatically.*

* Available on FF frame only.





LIGHTWEIGHT

FROM 6.7 KG TRANSPORT WEIGHT

How do you build a ready to ride folding wheelchair at just 6.7 kg and still deliver outstanding rigidity? It starts with advanced materials and the engineering expertise that defines QUICKIE.

The difference is something you feel instantly. **As the lightest QUICKIE aluminum folding frame**, every lift, transfer, and adjustment takes less effort.

TRANSPORT PROTECTION

Additional protection on the exposed axle receiver and backrest tube helps prevent damage when transporting or storing the wheelchair without wheels.



AEROSPACE-GRADE MATERIAL

Built from 7000 series aluminum, the frame delivers exceptional strength at a remarkably low weight. Compared to the 6000 series used in many wheelchairs, it enables thinner, lighter tubing without compromising durability.

INCREASED FRAME STIFFNESS

A refined frame design - including cross brace, axle stem, and castor links - increases stiffness by 56% compared to Xenon². More of your energy goes directly into forward motion.

OPTIMIZED TUBE DESIGN

The patented tube profile distributes loads efficiently and reduces stress points, limiting flex under pressure. Its oval shape minimizes material use, lowers weight, and prevents clamp rotation.

The result is a precise, responsive, and highly efficient ride.

PRECISION CONTROL

The three way adjustable caster system fine tunes angle, tracking, and height for precise control and stable tracking. Less correction. Less fatigue.



CONFIGURABILITY

BECAUSE EVERY BODY IS DIFFERENT

A wheelchair performs best when the frame matches your body and movement. The right setup improves posture and propulsion efficiency.

That's why the QS7 X series offers **20 frame configurations**, more than comparable folding wheelchairs, while staying light and rigid.

- **5 frame versions:** Fixed Front, Hybrid, Swing Away as well as Front Fold and Pro option
- **Multiple frame angles** to suit positioning, stability and centre of gravity
- **2 different frame insets** to support leg positioning and reduce front width
- **Adjustable backrest angle** for ergonomic posture (optional)
- **Rear seat height adjustment** to optimise propulsion and reduce strain
- **Centre of gravity adjustment** for precise, adaptable positioning

The **RE-DESIGNED CROSS BRACE** fits seamlessly beneath the seat and is barely noticeable. This creates a clean, minimalistic, and open-frame aesthetic you would normally only expect from a rigid wheelchair.





TAILORED CASTOR GEOMETRY

The castor -rear wheel relationship is key to efficiency. This is why each setup is optimized to your frame angle and inset for smoother handling and efficient propulsion.



QS7 X - THE MINIMALIST

Our lightest adjustable aluminium folding wheelchair. Sleek, minimalist design. The open fixed front frame simplifies car transfers. Fully adjustable for a personalised setup.



QS7 X HYBRID - THE STRONGEST

The reinforced fixed front frame delivers exceptional rigidity and supports users up to 140 kg. Its strength makes it ideal for use with power assist devices, while steeply angled tubes allow independent car transfers.



QS7 X-SERIES

A SERIES OF HIGH-PERFORMING FOLDING WHEELCHAIRS

There is a QS7 X just for you. Whichever of the four feature-rich models you choose, QS7 X-series brings a level of sophistication and elegance that ensures you'll always stand out.

QS7 X WITH FRONT FOLD - THE SMALLEST

Designed to integrate seamlessly into the front frame architecture.

QS7 X SA – THE ALLROUNDER

Features footrests that move to the side, making transfers more comfortable, especially from standing. Removing the footrests lets you position the wheelchair closer for easier transfers and reduces the overall pack size.



QS7 X PRO – THE EXPERT

The Pro frame features a fixed welded backrest and welded axle stems, reducing moving parts, cutting additional 300 g and maximises rigidity. Recommended for experienced riders who do not require later adjustment.





SMART OPTIONS

MAKE IT PERSONAL

NEW CARBON SIDEGUARDS

Longer shaped carbon sideguards, for better positioning and protection.



LED LIGHTS

Integrated LED lights blend seamlessly into the QS7 X design.

Three modes, high, low, and blink, improve visibility so you can see and be seen.
Easily detachable for recharging.

NEW LIGHTWEIGHT AXLE

Newly designed titanium lightweight axle saves 90 g compared to the standard axle.



NEW ALUMINIUM FOOTREST

Lighter, easier to assemble, and simple to adjust. Available as platform or divided, with easy angle, width, and depth adjustment.



OMOBIC

New compact aluminium brake from Omobic with smooth, easy operation, also suitable for reduced hand strength. The ergonomic grip and ribbed surface improve handling and stability.



SEATING

If you have specific seating requirements, you can easily upgrade your QS7 X with different JAY cushions and backrest. Did you know that you can also get your JAY J3 backrest in frame color? Your QS7 X will look even more stunning!

 **JAY**[®]


MODERN DESIGN

MADE PERSONAL

With its integrated caster forks, the “hidden” cross brace, and minimal visible hardware, **the QS7 X looks sleek and modern** from every angle.

Now it is up to you to give it your own personal touch:

- Choice of 22 colors for both frame and handrim
- Upholstery Trim Colors
- Anodized or powder coated castor fork options
- Anodized or powder coated wheel options
- Frame Graphic

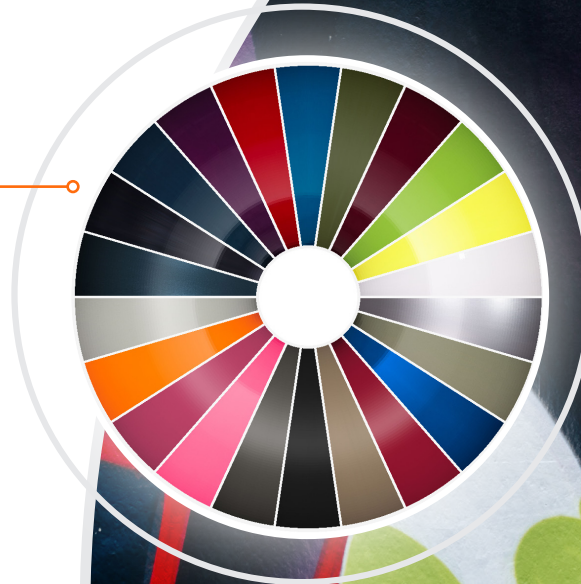

Empulse

THEY ARE GREAT TOGETHER:

The QS7 X is the perfect companion for any power assist device.



CREATE YOUR OWN **QS7 X**

WITH THE 3D VISUALIZER

Want to see what the QS7 X has to offer, how your QS7 X will look, and what it will weigh? Scan the QR code and build your QS7 X of your dreams in the **new 3D Visualizer**.



Scan the QR code
to learn more.

The QS7 X offers an extensive range of clinical-critical but also lifestyle options to suit a wide variety of needs.

For a complete overview of all available configurations, explore our new **Options & Accessories catalogue**.



Scan the QR code
to learn more.

Now you've read all about QUICKIE QS7 X, we think it's about time for you to experience it for yourself.

Book your test drive today!



built4me

Want something unique? Our Built-4-Me team can help you create a truly individual QS7 X that's out of this world! Ask your dealer for more information.

TECHNICAL DATA	QS7 X
Max. user weight	QS7 X FF - 110 kg, QS7 X Hybrid - 140 kg (136 kg with crash test option), QS7 X SA - 125 kg
Seat width	320 - 500 mm (in 20 mm steps)
Seat depth	340 - 500 mm
Front seat height	400 - 570 mm (SA) / 400 - 560 mm (FF)
Rear seat height	370 - 500 mm
Backrest height	250 - 475 mm
Backrest angle	from -15 to +13 / 75°-103°
Total width	SW +180/240/270 mm at 0°, 2°, 4° camber and 25" wheels The total width can exceed 800 mm in certain configurations.
Total length	1150 mm
Total height	1050 mm
Max. pivot width	760 mm
Max. safe slope for use of parking brakes	7°
Frame versions	Fixed Front, Swing Away, Hybrid, Front Fold, Pro
Frame angle	SA 70°/80°/90° (110°/100°/90° EU), FF 80°/88°/92° (100°/92°/88° EU)
Front wheel	3", 4", 5", 6", 7"
Rear wheel	24", 25"
Wheelchair weight	Swing Away: from 10,4 kg, Fixed Front: from 9,3 kg
Dismanteled weight	Swing Away: from 7,8 kg, Fixed Front: from 6,7 kg
Mass of heaviest part	2,1 kg
Removable parts	Rear wheels (for transportation)
Frame colours	22 colours
Crash tested	Yes (approved for transport in a vehicle according ISO 7176-19:2022)
Transport without user	Fix the wheelchair firmly and securely in the vehicle (e.g. car or plane)
Intended use and environment	For people unable to walk and with limited mobility; inside and outside. Operation by user or assistant.



Sunrise Medical GmbH

Kahlbachring 2-4

69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: international@sunrisemedical.de

For further information on the full specification, options, and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.



If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.eu.



Rev. A © 04.2026



www.SunriseMedical.eu

