




ZIPPIE VOYAGE

ADVENTURES
AWAIT

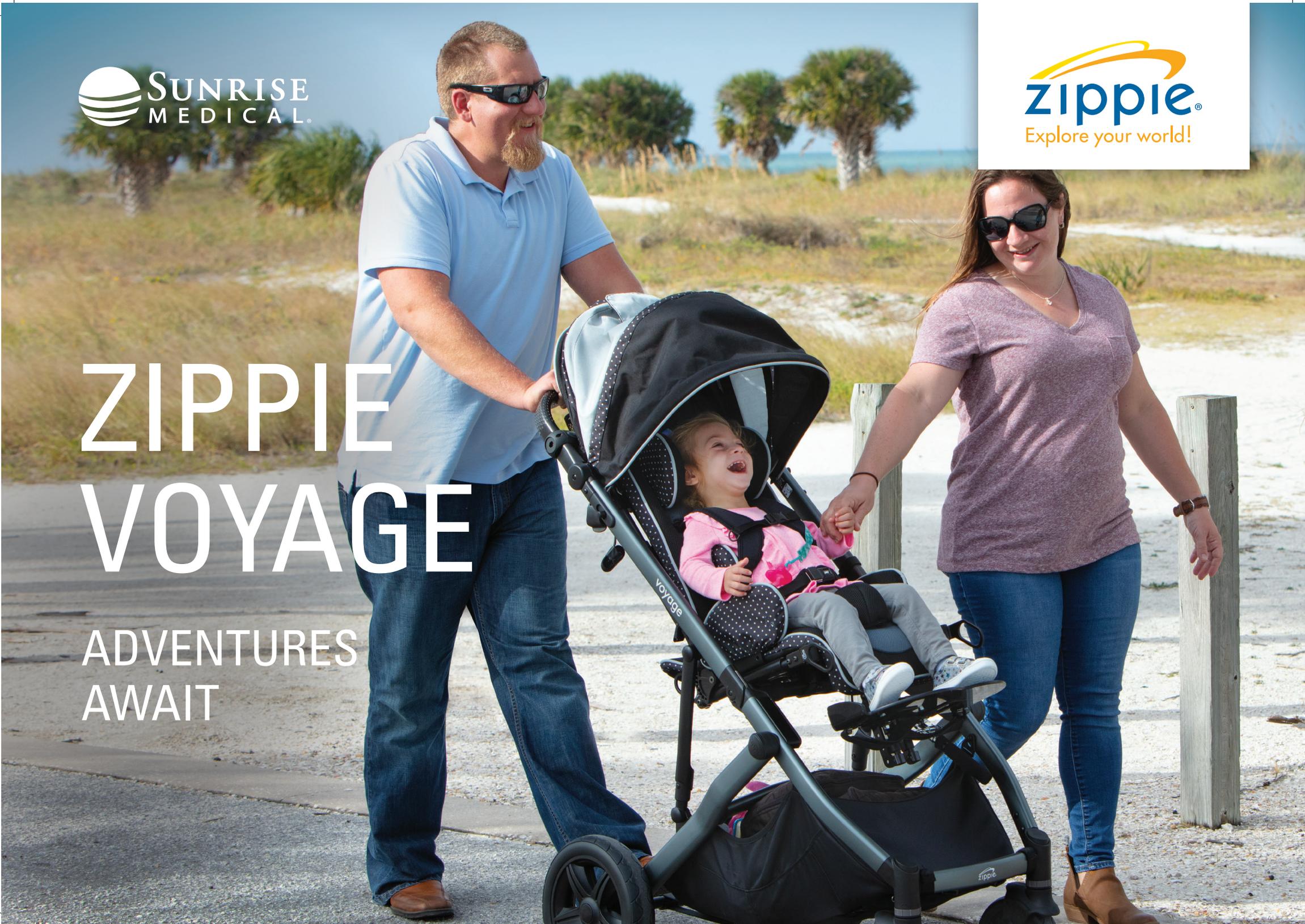

MODERN. EASY TO USE. BEST FOR BABY.



The universal Zippie Voyage early intervention device combines the convenience of a quick-folding lightweight stroller frame, Zippie's incredibly versatile seating for growing babies, and one of the largest selections of functional accessories in the industry! The result is an intuitive, easy-to-use stroller that makes no compromises when it comes to what's best for baby.

// CANOPY WITH MESH WINDOW
Keep baby shaded during the day and open the mesh window to check in.

// CARRY-ALL STORAGE BASKET
Basket conveniently store up to 7 kg of your essentials.

// SWIVEL FRONT WHEELS WITH SUSPENSION AND LOCK-OUT
Glide over and around obstacles while keeping baby calm, secure, and comfortable. Lock-out the front wheels for smooth tracking during long-distance strolls.



// ANGLE ADJUSTABLE STROLLER HANDLE
Responsive, one-hand steering that adjusts to the caregiver's height.

// MODERN STYLING
The Zippie Voyage base with its sleek styling and lightweight design has been specially designed to incorporate the Zippie Voyage seat system.

// FOOT OPERATED WHEEL LOCKS
Hands free! Lock wheels in place easily with your foot.

// AIRLESS WHEELS
Flat tires are a thing of the past with airless wheels.





// FORWARD (ANTERIOR) TILT

Assists with positioning to encourage achievement of developmental milestones such as bringing hands to the mouth, reaching for objects, or engaging in activities.

// REARWARD (POSTERIOR) TILT & RECLINE

Improves positioning to assist with feeding, digestion, respiratory function, and visual orientation.

**// HEIGHT & ANGLE
ADJUSTABLE FOOTREST**

Grows with your child and adjusts for various activities and positioning requirements throughout the day.



// SIMPLE TILT AND RECLINE ADJUSTMENT

Tilt and recline improves positioning and helps with feeding, digestion, respiratory function and visual orientation. Tilt (range 35 °) and recline (up to 180 °) the Zippie Voyage by simply using the easy-to-reach trigger handles.

SEATING FOR BABIES ON THE GROW!

When it comes to seating, the Zippie Voyage's number one priority is positioning your child for healthy growth and development from babyhood through their toddler years.



// HEAD SUPPORT
Height and width adjustable

// JAY CHEST STRAPS

// QUICK ADJUST TILT
Easily tilt forward or rearward for positioning needs.

// BACKREST
Height adjustable

// LATERAL TRUNK SUPPORT
Height and width adjustable

// HIP GUIDES
Height and width adjustable

// JAY PELVIC BELTS

// POMMEL
Available in two different versions

// FOOTREST
Height and angle adjustable footrest with optional foot straps and shoe holders

EASY HANDLING FOR DAILY USE

Quickly and comfortably attach and remove the seat from the stroller in just one simple motion. This makes transport easy and guarantees that the stroller fits in even a small car.



// INTUITIVE FOLDING

Easy-to-reach folding triggers fold the stroller in an instant into a handy pack size so you can quickly be on your way.

// REVERSIBLE SEAT

Face the seat towards you to interact with your child or forward so they can explore their surroundings.

// SLIDE-N-LOCK SEAT ADAPTER

Remove and re-attach the seat is as simple with the Slide-n-Lock adapter. Rails guide the seat into position until an audible lock indicates that the seat has been safely attached.



LIFE IS COLOURFUL!

Complement the modern design and attractive polka-dot upholstery of your Zippie Voyage with a selection of four distinctive colours.



FUNCTIONAL ACCESSORIES OF THE ZIPPIE VOYAGE

The Zippie Voyage is available with options and accessories designed to make daily tasks a breeze, so you can delight in the little joys of caring for your family.



// JCM HI-LO BASE

This optional hi-low base allows your little one to engage in activities at the table or on the floor, whilst experiencing all of the seating and positioning benefits.



// CUP HOLDER

For a water bottle or a coffee-to-go - this cup holder uses a self-leveling mechanism to prevent unwanted spills when strolling.



// BELLY BAR

The belly bar gives your child a safe and comfortable area to hold on to whilst strolling.



// UPPER EXTREMITY SUPPORT TRAY

This support tray is available in two versions – in light grey or with clear poly material. Both versions are depth and height adjustable.



// HOOK UTILITY

Useful hook to hold multiple bags. Easy to attach with velcro strap that stays put and will not slide around.



// UTILITY UNDER FRAME TRAY

Versatile hard tray storage option for e.g. batteries or vents. Maximum load up to 23 kg.



// FRONT MOUNT UTILITY PLATFORM

Front mount storage platform to enable easy access to carried objects.



// SHOE HOLDER

The shoe holders helps to improve overall positioning by holding the foot securely.

TECHNICAL DATA

Zippie Voyage // Early Intervention Device

Max. user weight:	34 kg	Turning Circle:	1780 mm
Seat width:	240 – 300 mm, with insert upto 155 mm (excl. upholstery)	Max. safe slope:	9° (for use of parking brake)
Seat depth:	200 – 330 mm (excl. upholstery)	Rear wheel:	12" Pneumatic with Airless Inserts
Lower leg length:	155 – 310 mm	Front wheel:	8"
Backrest height:	410 – 610 mm (excl. upholstery)	Seat weight:	9,8 kg
Backrest angle range:	80° to 180°	Base weight:	12,7 kg
Seat angle range:	35° (-5° to 30°)	Removable parts:	Seat and wheels - for transport only
Seat height:	575 – 640 mm (3 steps possible 575 / 600 / 640 mm)	Colour options:	4 colours - blue, red, beige or gray trims
Total width:	615 mm	Frame material:	Aluminium
Total length:	1100 mm	Crash tested:	Yes; tested in accordance with Annex A of ISO 7176-19 (2008)
Total height:	810 – 1220 mm	Transport without user:	Fixe the product safe in the vehicle (e.g. car or plane)
Folded length:	800 mm	Intended use and environment:	For children with limited mobility: indoor and outdoor use, operated by caregiver
Folded height:	310 mm		

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.eu or alternatively is available on request in large text.



Sunrise Medical GmbH

Kahlbachring 2-4

69254 Malsch/HD

Tel.: +49 (0) 72 53 / 9 80-0

Fax: +49 (0) 72 53 / 9 80-222

E-Mail: international@sunrisemedical.de

www.SunriseMedical.eu