



WHERE YOU WANT TO BE.

ELITE² FAMILY[®] SERIES

COMFORT, RELIABILITY
AND SAFETY ON THE ROAD



ELITE² Mini



ELITE² XS



ELITE² Plus



THE ELITE² FAMILY


An established scooter family with enhanced quality and upgraded service information - for improved lifetime and increased reliability. Easy to ride, affordable and environmentally friendly - choose the scooter that suits your needs.

The Elite² family combines many advantages so that you can feel safe and comfortable at all times. With the familiar ergonomic wig-wag control levers the Elite² family offers the much acclaimed hand / thumb control as a standard. On the LCD display you can easily read out the day's distance, total distance and driving speed. The correct sitting position can be achieved by adjusting the seat and steering column.

The Elite² scooters are the ideal travel companion with a maximum speed of 12 or 15 km/h (depending on model), which will bring you safely to your destination.

OPTIMAL SAFETY

- Integrated CE-marked clear and bright LED lights
- A secondary brake with electronic stop function
- Easily visible break lever with electronic stop function ensures instant braking if required
- Strong magnetic parking brake

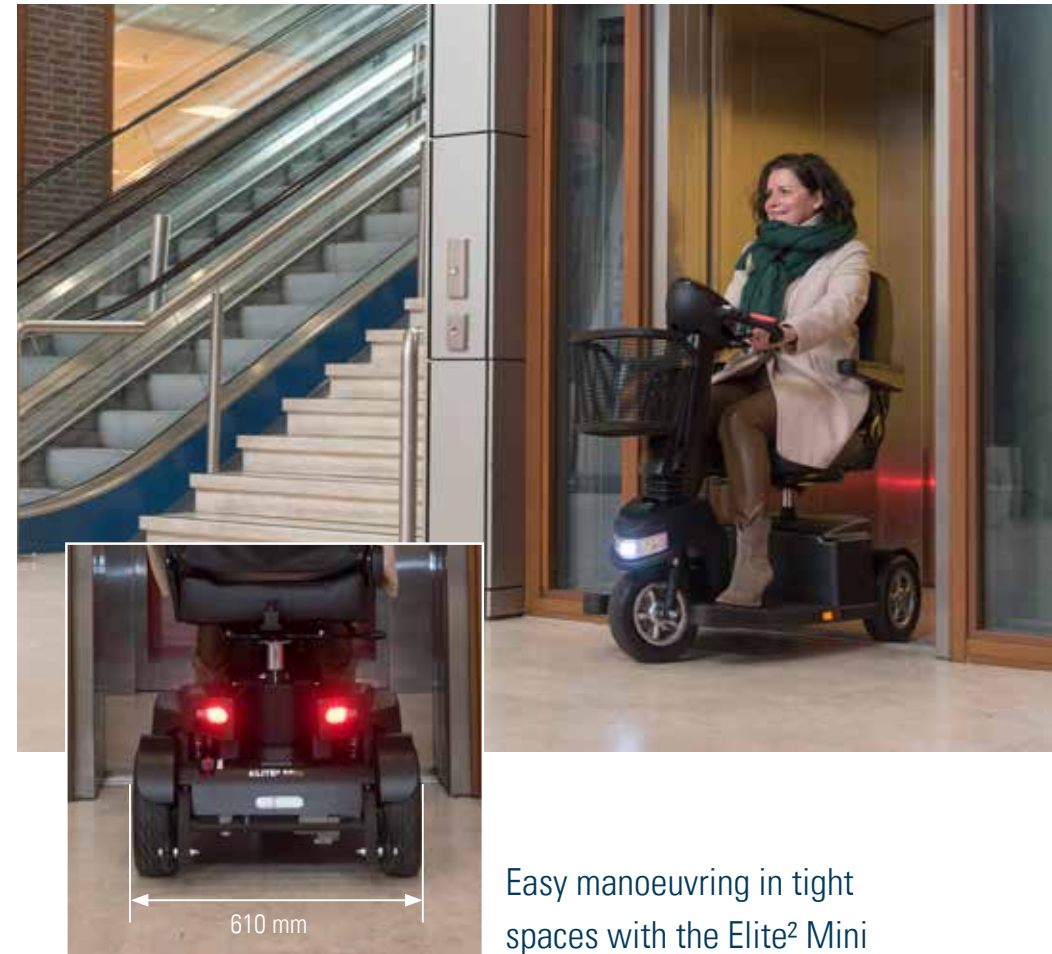

Easy manoeuvring in tight spaces with the Elite² Mini

ELITE² *Mini*



- Compact scooter for better manoeuvrability in narrow spaces/environments, in shopping malls, lifts
- Adjustable rear suspension for a comfortable drive for occasional use on a variety of surfaces

ELITE² *XS*



- Perfect balances ration between length and wheel size for easy driving
- Very comfortable drive thanks to bigger wheels and adjustable rear suspension for every user on a variety of surfaces

ELITE² *Plus*



- Extra comfortable drive thanks to biggest wheels and extended adjustable rear suspension
- Steers just like a car: light and easy to use
- Optional twist controller (hand throttle) available for more ergonomic hand position

HYBRID TYRES


Low risk of punctures due to hybrid tyres. The Elite² family series is characterised by the particularly robust hybrid tyres. This not only reduces costs for the user. The dealer can also profit enormously from this. Hybrid tyres are more stable and thus reduce the need for service appointments.

The tyres run much more smoothly and, especially on uneven ground, much more quietly as they can absorb unevenness well. It is therefore puncture-proof and makes the day trip easier - this increases safety and reduces the risk of getting stuck with a flat tyre.

A tyre technology with greater safety and comfort for carefree driving fun.

Available on Elite² XS and Elite² Plus.

Solid insert
Inner air tube

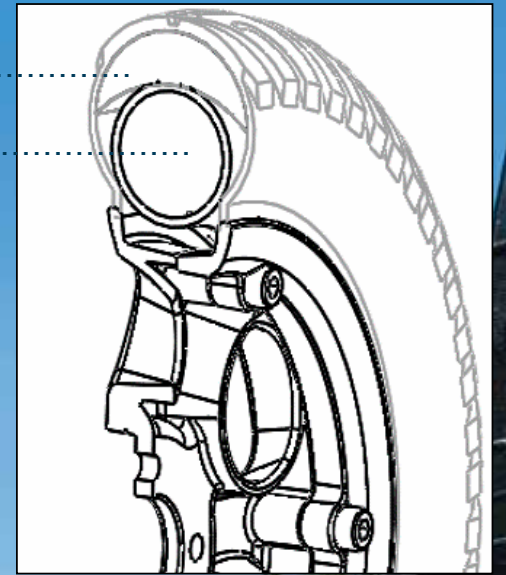

A close-up, low-angle shot of the front wheels of two scooters parked on a cobblestone path. The wheels are black with a deep, wavy tread pattern. The scooters' fenders and front forks are visible. The background is a blurred green lawn. The text is centered in the upper half of the image.

There is a wide range of options
and accessories that will provide
the scooter with all the
comforts, needs and your personal
preferences.

DETAILED FEATURES



Front bumper as standard to reduce damage (for Elite² Mini)

Brake lever in red for better and faster recognition

All Elite² Mini Scooters feature ergonomic soft touch keys

Full front bumper as standard to reduce damage to the front (e.g. against doors)

Hybrid tyres
(for Elite² Plus & Elite² XS)



EVERYTHING UNDER CONTROL



All Elite² Mini Scooters feature ergonomic hand thumb control.

The Elite² family scooters have ergonomically designed handles and the intuitive ergonomic finger and thumb operated wig-wag control. The handlebar is angle adjustable in to ensure comfort and good posture while driving.


The clear and easy to operate display is fitted with a reliable battery indicator and speedometer. Optional available an energy efficient LCD display with additional odometer and hour meter.



LUXURY SEATING

**The multifunctional options are
suitable for every trip**

- Easily adjustable back angle
- Height adjustable armrests
- Fully adjustable lumbar support
(on Elite2 XS and Elite2 Plus)
- Sliding and turning mechanism
for easy transfers

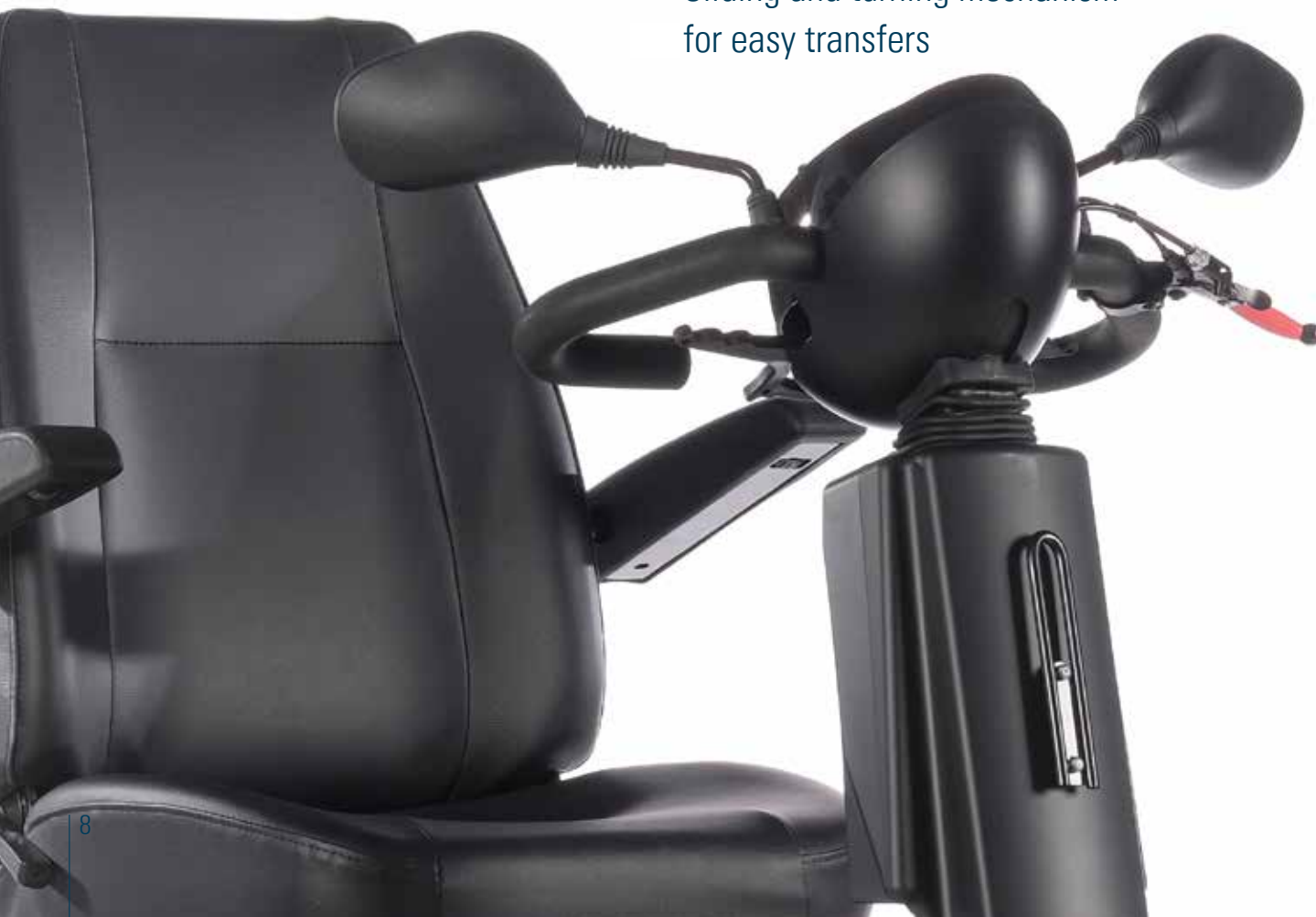

COMFORT ON EVERY TRIP

The Elite² family advanced seating solution takes comfort to the highest level – from the seat height, depth and recline adjustments, to the flip up, width, angle and depth adjustable comfort armrests. The armrest can be flipped back to facilitate transfers and can be also easily angle adjusted for optimum positioning. Every feature has been designed for comfort.



- Moving the cockpit



- 360° swivel seat enable comfortable transfers



- Depth adjustable seat



- Flip back armrests



Elite² Plus

ELITE² FAMILY[®] SERIES



DESIGN AND MAXIMUM SAFETY



Maximum safety on the road

To ensure the visibility on the road, the Elite² family scooters incorporates front, rear, indicators and brake LED lights. It also features mudguards, hand brake, anti-tip wheels and rear bumper as standard. All in a sturdy structure, approved for users up to 136 kg / 21 stone (Elite² Mini) and 175 kg / 27 stone (Elite² XS and Elite² Plus).



ELITE² FAMILY[®] SERIES



Comfort, reliability and safety on the **road**.

True to this theme, we have developed the Elite² family. It focuses on the improvement of existing scooters. The focus here is on rental and re-use. This enables dealers to keep existing parts - sustainable, sensible and cost-effective.



TECHNICAL DATA

	ELITE ² MINI	ELITE ² XS	ELITE ² PLUS
Type	3 wheel	3 wheel	3 wheel
Colour	Textured Black	Textured Black	Textured Black
Total Length (max.)	1270 mm / 50 inch	1390 mm / 54.7 inch	1390 mm / 54.7 inch
Total Width	610 mm / 24 inch	670 mm / 26.3 inch	650 mm / 25.6 inch
Seat Depth	440 mm / 17.3 inch	460 mm / 18.1 inch	460 mm / 18.1 inch
Seat Width	470 mm / 18.5 inch	475 mm / 18.7 inch	475 mm / 18.7 inch
Seat Surface Height	430 - 510 mm / 16.9 - 20.1 inch	440 - 540 mm / 17.3 - 21.2 inch	440 - 540 mm / 17.3 - 21.2 inch
Backrest Height	500 mm / 19.7 inch	510 mm / 20.1 inch	510 mm / 20.1 inch
Backrest Adjustment	Yes	Yes	Yes
Seat Adjustment (horizontal)	Yes	Yes	Yes
Maximum speed	10 km/h / 6.2 mph (via programming) / 12 km/h / 7.4 mph	15 km/h / 9.3 mph	15 km/h / 9.3 mph
Maximum Range	38,8 km / 24 miles	44 km / 27 miles	44 km / 27 miles
Maximum Safe Slope	9° at 136 kg / 21 stone	10° at 175 kg / 27 stone	10° at 175 kg / 27 stone
Turning Circle	2460 mm / 96.85 inch	2800 mm / 110.2 inch	2800 mm / 110.2 inch
Tyres	Pneumatic	Hybrid, Pneumatic or Solid	Hybrid, Pneumatic or Solid
Front wheel size	10.5"	13"	14"
Rear wheel size	10.5"	13"	14"
Total weight excluding batteries	70 kg / 11 stone	86 kg / 13 stone	94 kg / 14 stone
Total weight including batteries	110 kg / 17 stone	132 kg / 20 stone	140 kg / 22 stone
Battery type	40/60 Ampere	60 Ampere	60 Ampere
Maximum User Weight	136 kg / 21 stone	175 kg / 27 stone	175 kg / 27 stone



All over Europe, thousands of people have made their lives easier with a STERLING mobility scooter by Sunrise Medical. For over 30 years, STERLING have designed and built mobility scooters to the highest standards, giving the comfort and safety every customer needs. For further information on product options please check the order form, user manual and service manual on our website. You can view this document also in a larger font on a PDF-file at www.SunriseMedical.eu

Sunrise Medical GmbH

Kahlbachring 2-4
 69254 Malsch/HD
 Tel.: +49 (0) 72 53 / 9 80-0
 Fax: +49 (0) 72 53 / 9 80-222
 E-Mail: international@sunrisemedical.de
www.SunriseMedical.eu